



KATSAN

MEDICAL DEVICES

"Professional tools
for professional hands"

PRODUCT CATALOGUE INCLUDING;

SURGICAL SUTURES · ALCASORB ALCASORB RAPID
ALCALACTINE ALCALACTINE RAPID ALCALACTONE
ALCADINONE ALCAMIDE ALCALENE ALCASILK ALCASAFE
POLYESTER UHMWPE ANCHOR SUTURES SCREW
ENDOBUTTON ALCATER TAPE STERNUM STEEL

LAPAROSCOPIC PRODUCTS · LAPAFIX TROCARS PTIC
TROCARS VERESS NEEDLE SUCTION IRRIGATION SET
ENDOBAG

SURGICAL SUPPORT · PP MESH COMPOSIT MESH
HEMOSTAT SPONGE SKIN STAPLER BONE WAX

VETERINARY PRODUCTS CASSETTE SUTURES
ANTIBACTERIAL SUTURES





Katsan has a wide range of medical devices, that are being used in both human and veterinary practices. With over 40 years of experience, Katsan has been led to the highest priority in the field of quality control, the development of medical products and surgical sutures made in accordance with innovative technologies.

Founded in 1976, Katsan Medical Devices is one of the leading synthetic surgical suture manufacturers in Turkey. In addition to surgical thread production, Katsan has expanded its service and product range and also offers laparoscopic surgery supplies, and surgical supports devices as hemostats and surgical mesh.

Katsan has developed new production strategies with its investments in quality production and innovative technology; With the support of the project studies carried out by the Quality Assurance and R&D team, who are experts in their fields, in both human and veterinary fields, we provide CE and ISO 13485 certified high quality, sterile medical materials to many countries in North America, Europe, the Middle East, Asia, Africa and Oceania. It started to sell actively in the market and managed to quickly integrate into these markets with quick returns and creative solutions.



With our fully CE marked and ISO 13485 quality certificated product array that are produced in our T.R. Ministry of Health approved facilities, we produce world-class sterile surgical medical devices for a variety of surgical operations.

Our company success depends on high-performing and professional employees who are experts in their own fields and the creative solutions that we give to our customers.

At Katsan, we are responsible to improve quality of lives of the patients, meet the demands of customers and surgeons and we achieve these through ethical practice and by being respectful to our environment.

INCLUDES

WOUND CLOSURE PRODUCTS

Surgical Sutures • Barbed Suture
Orthopedic Sutures • Cervix Set • Sternum Steel

LAPAROSCOPY PRODUCTS

Lapafix Bladed & Bladeless
Trocars • Lapafix Optical Trocar Veress Needle • Suction
Irrigation Set • Endobag

SURGICAL SUPPORT PRODUCTS

PP Mesh • Composite Mesh • Hemostat
Sponge • Skin Stapler • Bone Wax

VETERINARIAN PRODUCTS

Cassette Sutures • Anti-Bacterial Sutures



WOUND CLOSURE PRODUCTS

Surgical Sutures • Barbed Sutures
Orthopedic Sutures • Cervix Set • Sternum Steel

SURGICAL SUTURES

Katsan Surgical Sutures, a leading synthetic surgical suture manufacturing company in Turkey, was founded in 1976 as a catgut producer. Katsan is striving for higher quality products and a bigger market share by using innovative technologies.

With our fully CE marked and ISO 13485 quality certificated product array that are produced in our T.R. Ministry of Health approved facilities, we produce world-class sterile surgical medical devices for a variety of surgical operations.



ABSORBABLE SUTURES



ALCASORB®
Polyglycolic Acid




ALCASORB® RAPID
Polyglycolic Acid Rapid



ALCALACTINE®
Poly (glycolic acid (90%)-co-lactic acid (10%))



ALCALACTINE® RAPID
Poly (glycolic acid (90%)-co-lactic acid (10%)) Rapid

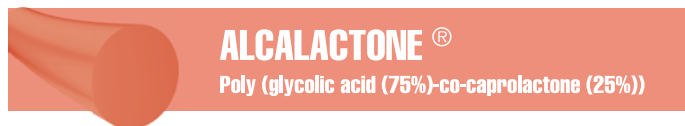

ALCALACTONE®
Poly (glycolic acid (75%)-co-caprolactone (25%))



ALCADINONE®
Polydioxanone

NON-ABSORBABLE SUTURES



ALCALENE®
Polypropylene




ALCAMIDE®
Polyamide



ALCASILK®
Silk


ANTIBACTERIAL SUTURES

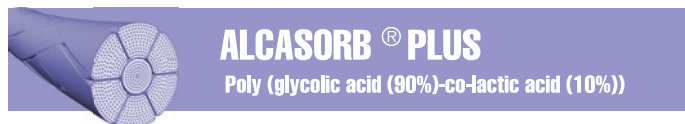

ALCASORB® PLUS
Poly (glycolic acid (90%)-co-lactic acid (10%))



ALCALACTINE® PLUS
Poly (glycolic acid (90%)-co-lactic acid (10%))

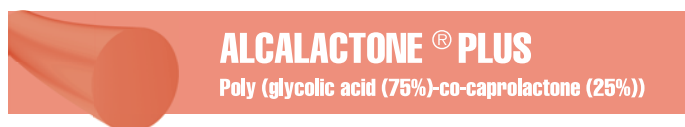

ALCALACTONE® PLUS
Poly (glycolic acid (75%)-co-caprolactone (25%))




ALCADINONE® PLUS
Polydioxanone

AREAS OF USAGE


Plastic surgery, urology, gastrointestinal surgery, skin closure, orthopedic surgery, cholecystectomy, hernia, laparoscopic incisions, peritoneum, fascia, ligation, general tissue approximation, emergency surgery, oral surgery, episiotomy, pediatric surgery, minor surgical procedures, subcutaneous, gastric surgery, intestine, colon, vesica urinaria, vaginal canal, uterus, etc.



ATRAUMATIC NEEDLE FORMS


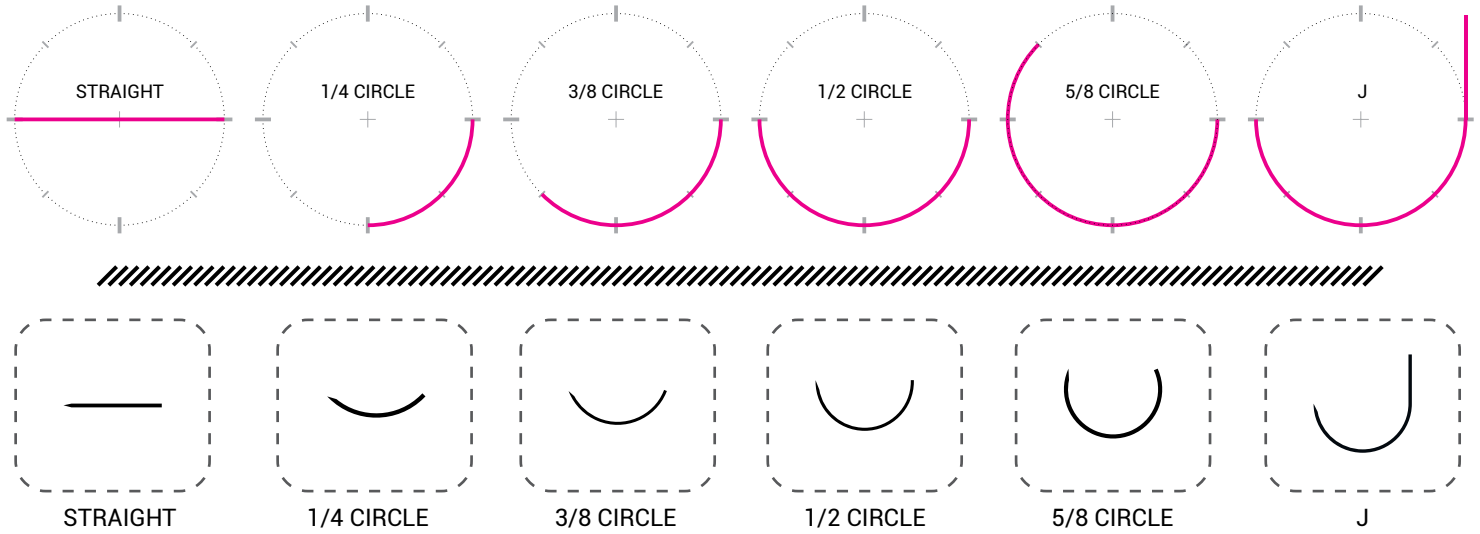


ATRAUMATIC NEEDLE TYPES

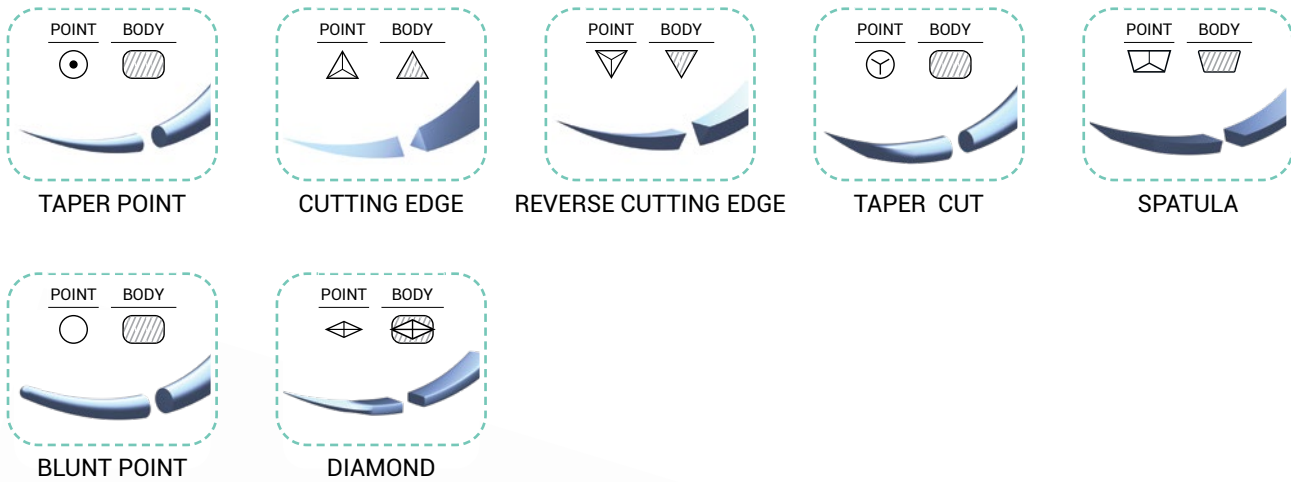

BARBED SUTURE

- Has microscopic etching which allows self-anchoring without the need of knots.
- The monofilament, synthetic material hydrolyzes within 180 days and tissue tension is maintained for about a month.
- Provides less blood loss during the surgery
- Less surgical difficulty
- Given its monofilament structure, reduces the local inflammation and risk of infection
- Allows the surgeon to complete the procedure easily and speeds up the process
- Minimal sacrifice of closure strength and scar cosmesis

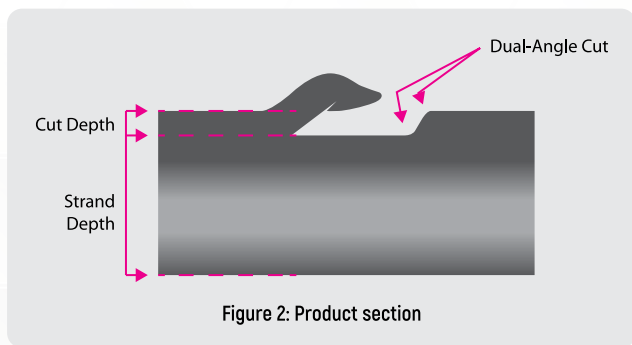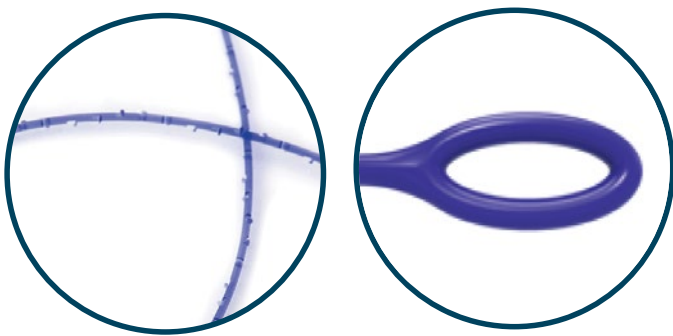

USP	Metric
USP 2	Metric 5
USP 1	Metric 4
USP 0	Metric 3,5
USP 2/0	Metric 3
USP 3/0	Metric 2
USP 4/0	Metric 1,5

Table 1: USP/METRIC Ranges





Product Specification	Technical Information
Color	Clear (undyed) and violet (dyed) colored monofilaments (D&C Violet No.2)
Chemical Structure	Polydioxanone (100 %) - synthetic, absorbable, monofilament surgical suture
Wound Support Profile	~60 days
Absorption Profile	~180 - 210 days
Content of residues	Monomers 1.0 w/w % Catalysts 100 ppm Heavy metals 10 ppm Solvents 500 ppm
Sterilization Method	EO (Ethylene Oxide)
Product Type	Uni-directional spiral type
Moisture Content	400 ppm
Retention ratio of tensile strength	About 50 % of the initial strength after two weeks
Storage Condition	Low temperature (1~30°C) Low humidity (less than 60%RH) Dry chamber or the refrigerator at -10°C

Chart 1: Technical Specifications

Size		Distance (mm)	
USP Size	EP (Gauge No.)	Minimum average value	Maximum average value
4-0	1.5	0.360	0.440
3-0	2	0.410	0.490
2-0	3	0.460	0.540
0(1-0)	3.5	0.510	0.590
1	4	0.560	0.640
2	5	0.610	0.690


Chart 2: Barbed Spacing (Distance between barbs according to USP/ METRIC range)

ORTHOPEDIC SUTURES

Orthopedic surgery is a medical specialty that focuses on treating activity-related orthopedic problems. These types of injuries usually involve muscles, bones, joints, and ligaments and often occur during athletes' intense physical activities. Sports surgeons use surgical and non-surgical methods to identify, treat and prevent these types of injuries.


UHMWPE SUTURES

ALCATER ST
Ultra high molecular weight polyethylene (UHMWPE)

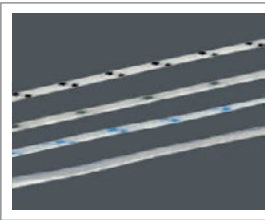

USP 2
Multiple size and colors is ideal for complex multi-strand procedures where better visibility and differentiation is a must.

ALCATER ST USP SIZE CHART							
COLORS	0/2	0	2	3&4	5	7	10
NATURAL WHITE	X	X	X	X	X	X	X
BLACK / WHITE	X	X	X	X	X		
BLUE / WHITE	X	X	X	X	X		
GREEN / WHITE	X	X	X	X	X		

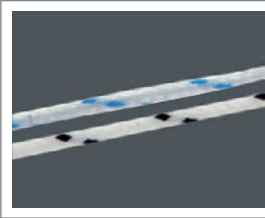

Approximation and fixation of soft tissue, cardiovascular surgery, use of allograft tissue for orthopedic procedures.



ALCATER ST TAPE
Ultra high molecular weight polyethylene (UHMWPE)

ALCATER ST TAPE USP SIZE CHART				
COLORS	1/4 mm	2 mm	3 mm	5 mm
NATURAL WHITE	X	X	X	X
BLACK / WHITE	X	X	X	
BLUE / WHITE	X	X	X	
GREEN / WHITE	X			

ALCATER ST TAPE
gently counters and lays flat against surface,


ALCATER ST TAPE will not pull through tissue as easily as sutures,



Ideal for applications where tissue pull through may be a concern,

Broad footprint for a uniform level of pressure over the surface area.



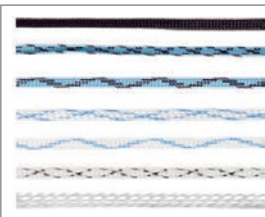


STRONGER THAN SUTURE INDICATION
The flat surface of the tape helps brace and distribute more consistent loads while still gaining the strength. Fixation of soft tissue into the bone. The use of allograft tissue for orthopedic procedures. Available in a variety of widths ranging from 1 mm to 5 mm.



Suture Anchor

Suture anchor intended to assist arthroscopic and orthopaedic reconstructive procedures that require soft tissue fixation due to injury or degenerative disease. Shoulder Indications; Bankart repair, SLAP lesion repair, Acromioclavicular (AC) separation, Rotator cuff repair, Capsule repair or capsulolabral reconstruction and Biceps tenodesis.

Advantages: Designed to securely re-attach soft tissue to bone. Available in two size for arthroscopic surgery: 3 mm and 5 mm. Minimize the possibility of tendon rupture. Open space inside of screw helps bone growth to allow bone marrow and blood flow.

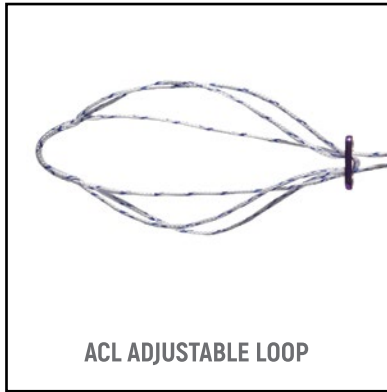

ACL Adjustable Loop

ACL Adjustable Loop is a self-latching adjustable button loop. It is made of ultra-high molecular weight polyethylene (UHMWPE) non-absorbable polymer. Its sizes change between 12 and 60 mm.

Advantages:

- Provides excellent durability and strength.
- According to the design of unparalleled 3-point locking, a suture is easily placed and the suture creep is decreased.
- One size fits all your ACL loop repair diameters.
- It has a high pull strength.
- The shiny titanium button provides a high visibility.
- To reduce the coiling of the loop and suture; Packaging is optimized by using a single, long, protective card.




ACL Closed Loop

ACL Closed Loop is a self-latching adjustable loop button for Femoral and Cortical Fixation.

Made of ultra-high molecular weight polyethylene (UHMWPE) non-absorbable polymer. Its sizes change between 12 and 60 mm

Advantages:

It has high tensile strength and durability.
The button part is shiny titanium.

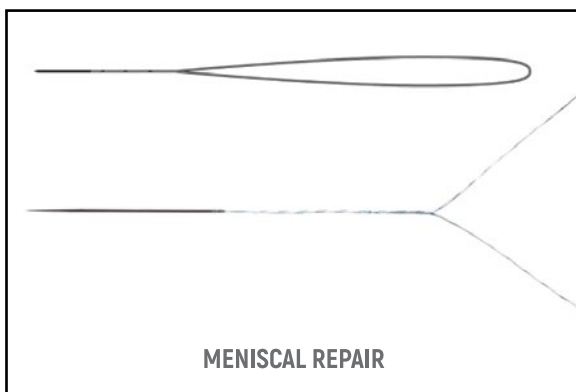

Meniscal Repair Sutures

Meniscal Repair Sutures which are made of an ultra-high molecular weight polyethylene (UHMWPE) non-absorbable suture. The sutures are commonly used in grafting, ACL repair, rotator cuff repair, achilles tendon repair, other soft tissue and anchor applications.

Advantages:

- Excellent for grafting, ACL repair and rotator cuff repair, achilles tendon repair, Other soft tissue and anchor applications.
- Suitable for femoral cortical fixation.
- Available in several variations for shuttling tendons: Closed loop, Split stitch, Shuttle loop, Passing loop.

ANCHOR SUTURES & SCREWS



Product	Suture	Needle
Titanium Suture Anchor Ø 3mm	90cm x 1	2 needle
Titanium Suture Anchor Ø 5mm	90cm x 1	2 needle
Titanium Suture Anchor Ø 3mm	90cm x 2	4 needle
Titanium Suture Anchor Ø 5mm	90cm x 2	4 needle

Suture color options: green, White, White/green, blue, White/blue
Suture type: UHMWPE



CERVIX SET


Cervix set is a specialized medical kit designed for use in gynecological procedures involving the cervix, which is the lower part of the uterus that opens into the vagina. The set typically includes various surgical instruments and devices needed to perform diagnostic and therapeutic procedures on the cervix. Cervix sets are essential for safe and effective gynecological care, ensuring that healthcare providers can perform procedures accurately and with minimal discomfort to the patient.



STERNUM STEEL

Sternum steel suture is a type of surgical suture made of stainless steel wire, specifically designed for use in closing the sternum (breastbone) after surgeries such as open-heart surgery. These sutures are durable, strong, and biocompatible, making them ideal for securing the sternum, which needs to be stable during the healing process after being split open during procedures like a median sternotomy. After the surgery, these steel

sutures are used to tightly close the sternum to ensure proper alignment and stabilization during the healing process. The strength of the steel wire helps keep the sternum securely in place while the bone heals, preventing complications such as sternal dehiscence (separation of the sternum).

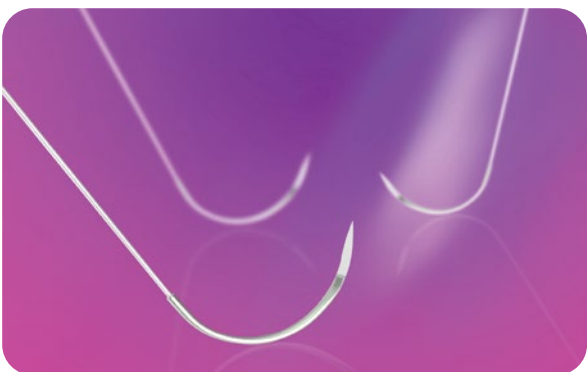

Suture	Suture Characteristics
Colour	Steel Grey
Raw Material	Stainless steel 316 LVM
Material Type	Monofilament
Absorption	Provides prolonged tensile strength in tissue & presence of infection
Characteristics	High tensile strength, good tissue compliance
Application	In general soft tissue approximation and/or ligation, including use in cardiovascular and orthopedic procedures.
Tensile Strength	High tensile strength.
Knotting Strength	Excellent
Tissue Reaction	Minimal
Sterilisation	Ethylene oxide
Shelf Life	5 years
USP Sizes	10/0 to 7
Packaging	12pcs / 24pcs or 36pcs in one box



LAPAROSCOPY PRODUCTS

Lapafix Bladed & Bladeless
Trocars • Lapafix Optical Trocar Veress Needle • Suction
Irrigation Set • Endobag


LAPAROSCOPIC TROCARS

LAPAFIX® trocars are a medical device which is placed in the body from a few different points without making wide incisions in the abdomen and chest region, one end is inside the abdomen, and the other end is located outside the abdomen, and can reach the patient's body without incision. The camera and the light source can be sent to the abdominal cavity through these trocars placed at several points, and the surgeon can follow the abdomen from the screen. Other hand tools and other necessary tools are transmitted and the operations are performed.

LAPAFIX® bladed trocar is available in ribbed fixation cannula configuration and standard trocar length (100mm) with 5 mm, 10 mm, 11 mm, 12 mm diameters.

LAPAFIX® bladed trocar has a sharp linear blade and blade shield with a spring-loaded locking mechanism. After penetration into the abdominal or thoracic cavity, the shield advances to cover the linear blade and reduces the possibility of injury in intra-abdominal or intra-thoracic structures. The trocar cannula contains an internal seal to prevent gas leakage when operating instruments are inserted or removed without loss of pneumoperitoneum. There is a stopcock valve for gas insufflation and desufflation. LAPAFIX® bladed trocar is sharper than reusable trocars and generally requires less force for insertion. Excessive force may reduce the user's control of the angle and the depth of insertion, which may increase the risk of injury to internal structures.

BLADED TROCARS

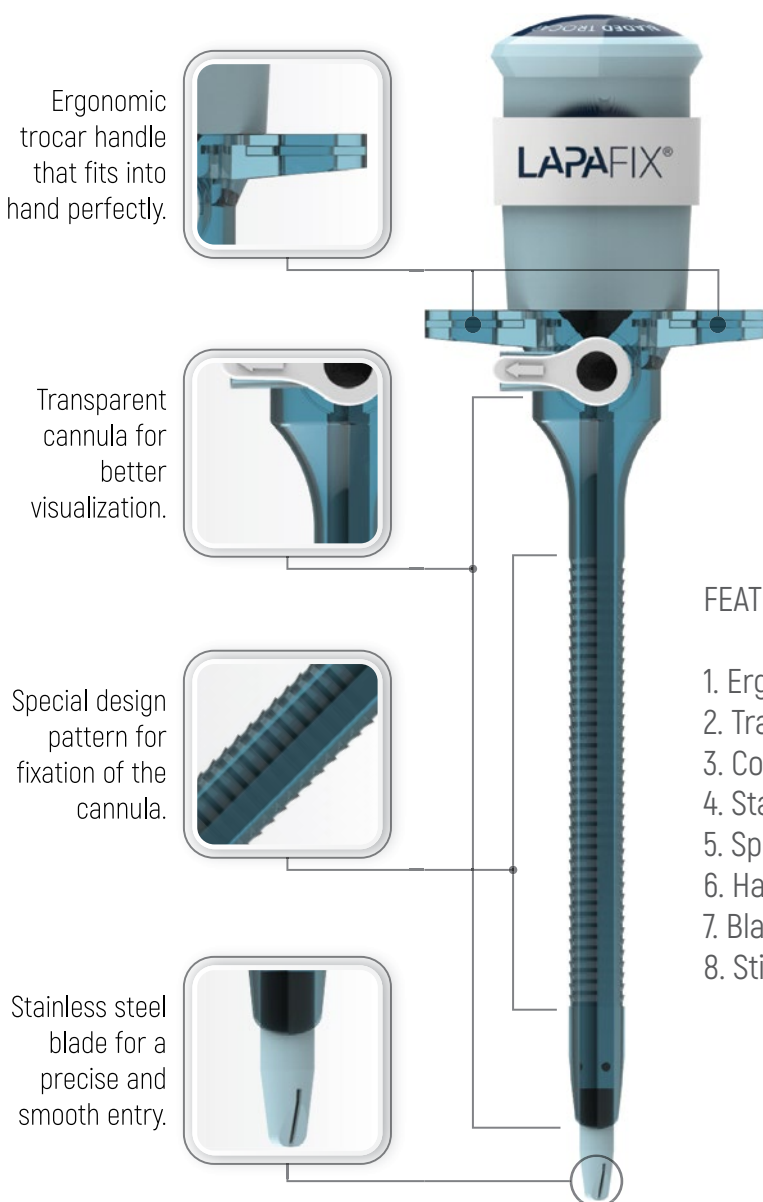


LAPAFIX®
CONFIDENT ACCESS, SAFE SURGERY

FEATURES

1. Ergonomic trocar handle that fits into hand perfectly.
2. Transparent cannula for better visualization.
3. Compatible with 5 mm instruments.
4. Stainless steel blade for a precise and smooth entry.
5. Special design pattern for fixation of the cannula.
6. Hand pressed dynamic locking for safety shield.
7. Blade position Indicator.
8. Stitching holes for fixation purposes.

BLADELESS TROCARS

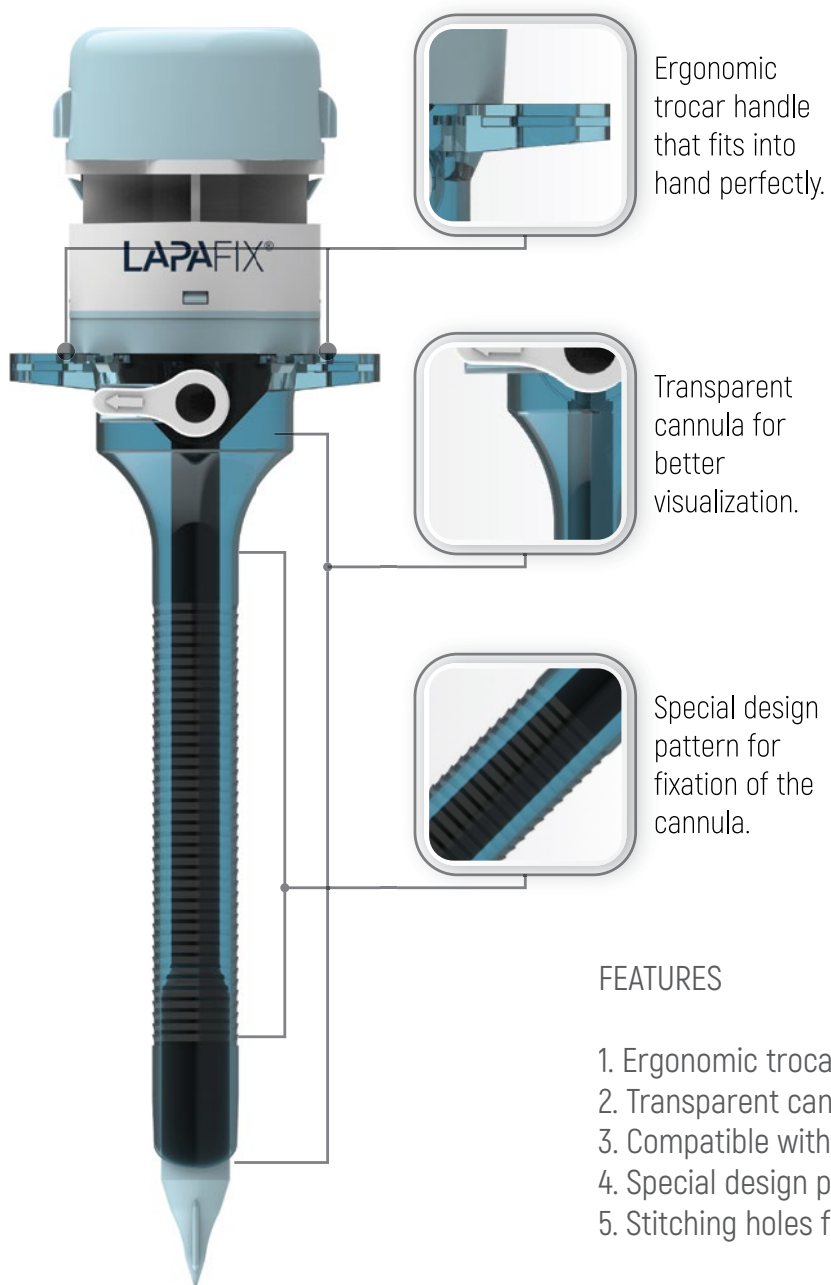

FEATURES


1. Ergonomic trocar handle that fits into hand perfectly.
2. Transparent cannula for better visualization.
3. Compatible with 5-10 mm instruments.
4. Special design pattern for fixation of the cannula.
5. Stitching holes for fixation purposes.



LAPAFIX®
CONFIDENT ACCESS, SAFE SURGERY



OPTICAL BLADELESS TROCAR

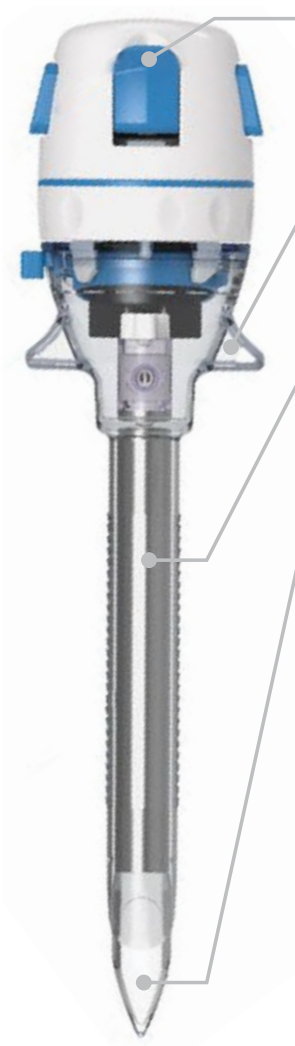

Laparoscope Locking Tab

- Holds camera locked in place during visually guided entry
- Helps provide visual clarity and protects camera from scratching or falling out

Suture Handle

- Fasten suture to the handle to fix the trocar in abdomen

Stability Thread

- Increases abdominal wall retention to avoid trocar slip-outs

Bladeless Tip

- Reduced need for fascial closure
- May result in a low port site herniation rate
- Separates, rather than cuts, along tissue fibers, pushing tissue and vessels away
- Reduced risk of abdominal wall vasculature injuries with bleeding complications

Universal Seal

- No extra reducer caps needed
- Maintains pneumoperitoneum
- Accommodates instruments sizes from 5mm to 12mm

Funnel Lead-In

- Sloped funnel enables instrument to insert smoothly

Stopcock Valve

- 90° rotation to charge or discharge CO2



Optical Entry

- Enables visualization of tissue layers during insertion

Cap Release Switch

- Release the cap to facilitate large specimen removal

PRODUCT NAME	ODEL	SPECIFICATION	DESCRIPTION
Optical Bladeless Trocar (Single Pack)	OT15M	15mm x 100mm	One obturator with one cannula
Optical Bladeless Trocar (Double Kit)	OT15MD	15mm x 100mm	One obturator with two cannulas


VERESS NEEDLE



Key Features of KATSAN Veress needle

- Transparent handle material allowing doctors to clearly see the needle indicator.
- Ribbed handle structure aiding easy penetration of the needle into the abdominal cavity.
- Ergonomic design for comfortable use.
- Appropriate channel diameter for quick CO2 loading.
- Unidirectional air valve preventing misuse.
- Spring-loaded blunt tip mechanism to protect the patient from injuries.
- Needle indicator with a distinctive red color for easy differentiation during penetration.

Product	Ref No	Features	Packaging
VERESS NEEDLE (2.1mm X 120mm)	EVN 12	120 mm	25 PCS / BOX
VERESS NEEDLE – LONG (2.1mm X 150mm)	EVN 15	150 mm	25 PCS / BOX

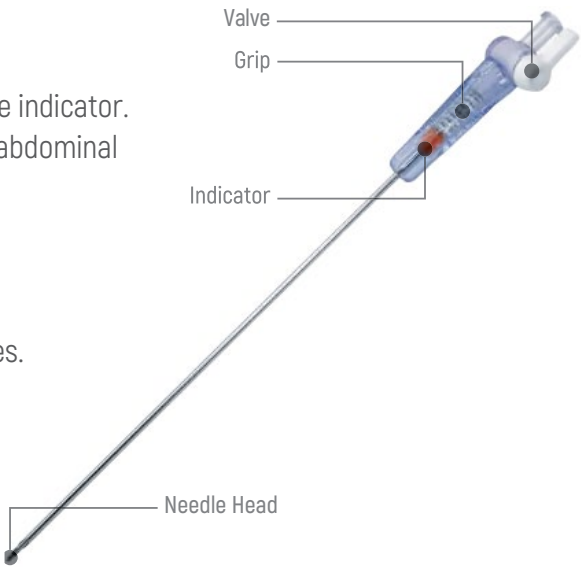

The Veress needle is a medical instrument used at the beginning of laparoscopic surgical procedures to inflate the abdominal area (pneumoperitoneum) and provide access to the abdominal cavity. The Veress needle is the first step in gaining access to the abdominal cavity. After the needle is inserted into the abdominal wall, the abdominal cavity is expanded by insufflating it with a gas such as air or carbon dioxide (CO2). This expansion allows the surgeon to better visualize internal organs and facilitates the more comfortable use of surgical instruments. While inserting the needle, the spring-loaded rounded tip (safety mechanism) moves inward, piercing through the outer layers of the abdomen. Once the needle enters the abdominal cavity, the rounded tip (safety mechanism) retracts, preventing injuries to internal organs.

SUCTION IRRIGATION SET

- Made of medical-grade PVC for safety and flexibility.
- Dual-entry design allowing the use of two fluid bags simultaneously.
- Differentiation of suction/irrigation channels with a blue line.
- Clips for regulating liquid flow.
- Enables the doctor to infuse medication when necessary.
- Provides the opportunity for specimen collection when needed.
- Ergonomic Design for Comfortable Use
- Different colors for suction/irrigation buttons for easy differentiation.
- Raised design to assist in distinguishing suction/irrigation buttons by touch.
- Spraying/Standard Modes to Accommodate Different Doctor Preferences
- Offers different modes to meet the preferences of various doctors.
- Constructed from high-density SUS304 material.
- Design prevents light reflections during operations.
- Perforated tip structure prevents clogging and accelerates irrigation.

Product	Ref No	Features	Packaging
Suction Irrigation Set	ESI05M	5mm X 330mm	10 PCS / BOX
Suction Irrigation Set	ESI05L	5mm X 450mm	10 PCS / BOX



STANDART ENDOBAG

Inner Shaft

- Push the endobag out of the introducer to make operation more convenient

Closure String release switch

- Opening the switch to detach the closure string from the introducer

Introducer

- Medical grade ABS material to make patients more safe
- 33cm length of introducer to reach deeper target areas

Memory Wire

- It could open automatically to simplify the operation and save doctor's time
- Adopting Ni-Ti alloy as material with high flexibility

Endobag

- Medical grade TPU/Nylon material are optional to meet different doctors'demands
- High intensity to ensure no damage of endobag
- High water resistance to ensure no pollution of liquid leakage
- Multi-volume for different sizes of specimen

Product	Model	Spesification	Packaging
Detachable Endobag with 100ml TPU material	EB100DT	10.2mmx230mm	10 pcs/box
Detachable Endobag with 200ml TPU material	EB200DT	10.2mmx230mm	10 pcs/box
Detachable Endobag with 400ml TPU material	EB400DT	10.2mmx230mm	10 pcs/box
Detachable Endobag with 1000ml TPU material	EB1000DT	12.2mmx230mm	10 pcs/box
Detachable Endobag with 200ml nylon material	EB200DN	10.2mmx230mm	10 pcs/box
Detachable Endobag with 400ml nylon material	EB400DN	12.2mmx230mm	10 pcs/box



SIMPLE ENDOBAG

Inner Shaft

- Adopting medical grade PVC as material to ensure safety and flexibility

Introducer

- Medical grade ABS material to make patients more safe
- 20cm length of introducer to reach any target areas

Memory Wire

It could open automatically to simplify the operation and save doctor's time

Adopting Ni-Ti alloy as material with high flexibility

Endobag

- Medical grade TPU/Nylon material are optional to meet different doctors'demands
- High intensity to ensure no damage of endobag
- High water resistance to ensure no pollution of liquid leakage
- Multi-volume for different sizes of specimen

Product	Model	Spesification	Packaging
Simple Endobag with 100ml TP bag material	EB100ST	10.2mmx230mm	10 pcs/box
Simple Endobag with 200ml TPU bag material	EB200ST	10.2mmx230mm	10 pcs/box
Simple Endobag with 400ml TPU bag material	EB400ST	10.2mmx230mm	10 pcs/box
Simple Endobag with 800ml TPU bag material	EB800ST	12.2mmx230mm	10 pcs/box
Simple Endobag with 1000ml TPU bag material	EB1000ST	12.2mmx230mm	10 pcs/box
Simple Endobag with 200ml nylon bag material	EB200SN	10.2mmx230mm	10 pcs/box
Simple Endobag with 400ml nylon bag material	EB400SN	12.2mmx230mm	10 pcs/box



SUPER SIMPLE ENDOBAG

Inner Shaft

- Medical grade ABS material to make patients more safe
- 20cm length of introducer to reach any target areas

Memory Wire

- It could open automatically to simplify the operation and save doctor's time
- Adopting Ni-Ti alloy as material with high flexibility

Endobag

Medical grade TPU/Nylon material are optional to meet different doctors'demands
High intensity to ensure no damage of endobag
High water resistance to ensure no pollution of liquid leakage
· Multi-volume for different sizes of specimen

Indicator

- Ergonomic design to avoid the string from falling into the abdomen through trocar

Product	Model	Spesification	Packaging
Simple Endobag with 100ml TP bag material	EB100UT	10.2mmx230mm	10 pcs/box
Simple Endobag with 200ml TPU bag material	EB200UT	10.2mmx230mm	10 pcs/box
Simple Endobag with 400ml TPU bag material	EB400UT	10.2mmx230mm	10 pcs/box
Simple Endobag with 800ml TPU bag material	EB800UT	12.2mmx230mm	10 pcs/box
Simple Endobag with 1000ml TPU bag material	EB1000UT	12.2mmx230mm	10 pcs/box
Simple Endobag with 200ml nylon bag material	EB200UN	10.2mmx230mm	10 pcs/box
Simple Endobag with 400ml nylon bag material	EB400UN	12.2mmx230mm	10 pcs/box





SURGICAL SUPPORT PRODUCTS

PP Mesh • Composite Mesh • Hemostat
Sponge • Skin Stapler • Bone Wax


MESHES

Meshes are usually made from durable and biocompatible materials such as polypropylene, polyester, or polyethylene. They are utilized during surgery to repair weakened or torn tissues or to provide support to areas affected by disease. For instance, in hernia surgeries, meshes can be used to strengthen the abdominal wall or repair a torn muscle.

Meshes are designed to support the tissue healing process and reduce potential complications. They encourage tissue growth while providing a resilient structure to facilitate healing.

PP Mesh

Teramesh® Non-Absorbable

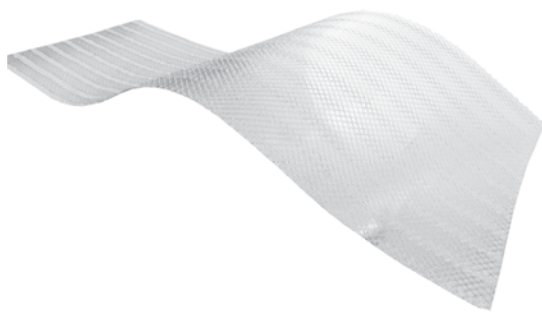

TERAMESH® Non-Absorbable is made of polypropylene non-absorbable monofilament fibers. This mesh implant is used for abdominal wall stabilization in cases such as hernia. After implantation the PP mesh remains stabilized and provides ease of usage as far as is not degraded in the body.

TERAMESH® Non-Absorbable is sterilised by ethylene oxide.


TERAMESH® Non-Absorbable is biocompatible and has ISO 10933 biocompatibility test reports.

PRODUCT	DIMENSIO	REF CODE	PACKAGING
TERAMESH® NON-ABSORBABLE	6cm x 11cm	FM0611	10 Pcs / Envelope
TERAMESH® NON-ABSORBABLE	10cm x 12cm	FM1012	10 Pcs / Envelope
TERAMESH® NON-ABSORBABLE	10cm x 15cm	FM1015	10 Pcs / Envelope
TERAMESH® NON-ABSORBABLE	15cm x 15cm	FM1515	5 Pcs / Envelope

PRODUCT	DIMENSIO	REF CODE	PACKAGING
TERAMESH® NON-ABSORBABLE	15cm x 30cm F	M1530	5 Pcs / Envelope
TERAMESH® NON-ABSORBABLE	30cm x 30cm F	M3030	5 Pcs / Envelope
TERAMESH® NON-ABSORBABLE	7,5cm x 15cm F	M7515	10 Pcs / Envelope

COMPOSITE MESH

Teramesh® Composite

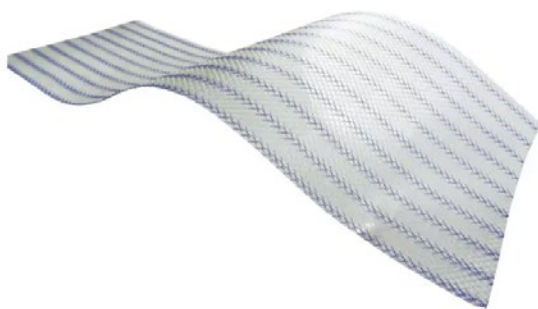

Sufficient Strength

Before absorption – rigid enough for easy manipulation
After partial absorption – retains sufficient strength for enduring normal abdominal pressure.

Less Stiff

After partial absorption, lower risk of pain with more soft and pliable mesh.

TERAMESH® Composite is a partially absorbable knitted mesh. Manufactured from approximately equal parts of absorbable Polyglycolic acid – co – caprolactone – 25 monofilament fiber and non-absorbable polypropylene monofilament fiber. Absorbable part of mesh is homogeneously absorbed in 90 – 110 days, after absorption of PGCL component, only PP mesh remains and provides maximum patient comfort.

TERAMESH® Composite can be used in open surgery and laparoscopic operations.

TERAMESH® Composite is sterilised by ethylene oxide.

TERAMESH® Composite is biocompatible and has ISO 10933 biocompatibility test reports.

Minimized foreign body

Lower risk of long – term side effects due to less residual foreign body remained after partial absorbability.

Its wide pore structure ensures maximum integration with the tissue and increases biocompatibility and ease of use

It has purple colored lines that will increase visibility in use.

PRODUCT	DIMENSIO	REF CODE	PACKAGING
TERAMESH® COMPOSITE	6cm x 11cm	PAM0611	10 Pcs / Envelope
TERAMESH® COMPOSITE	10cm x 12cm	PAM1012	10 Pcs / Envelope
TERAMESH® COMPOSITE	10cm x 15cm	PAM1015	10 Pcs / Envelope
TERAMESH® COMPOSITE	15cm x 15cm	PAM1515	5 Pcs / Envelope

PRODUCT	DIMENSIO	REF CODE	PACKAGING
TERAMESH® COMPOSITE	15cm x 30cm	PAM1530	5 Pcs / Envelope
TERAMESH® COMPOSITE	30cm x 30cm	PAM3030	5 Pcs / Envelope
TERAMESH® COMPOSITE	7,5cm x 15cm	PAM7515	10 Pcs / Envelope

HEMOSTATS



Hemostat is a medical tool used to stop or control bleeding. Hemostats are typically used during surgical procedures or in the treatment of traumatic injuries to control blood flow. These tools work by compressing blood vessels or other structures causing the bleeding.

TERACEL® is an absorbable and sterile woven fabric created by the oxidation of regenerated cellulose. The agent in question is light yellow in color and has a caramel-like odor. It should be stored at controlled room temperature. There may be a slight change in color over time, but this does not affect the performance of **TERACEL®**.

TERACEL® Original Haemostat is an absorbable and sterile knitted fabric, which is generated by oxidation of regenerated cellulose. Said agent is white with a pale yellow color and caramellike smell.

TERACEL® Original Haemostat is used for controlling capillary bleedings and other kinds of bleedings during surgical operations. TERACEL® Original Haemostat can be suitable for use where other kinds of haemostatic applications are not appropriate or convenient.

- Excellent Haemostasis in 2 – 3 minutes
- Low pH has antimicrobial effect
- Very good handling characteristics
- Does not stick to instruments
- Typically dissolves in 7 – 14 days
- Bactericidal properties
- Compliance with the requirements of USA
- Plant based product with minimum tissue reaction
- Do not cause irritation to applied body/tissue
- TERACEL® Original Haemostat are sterilized by Gamma irradiation technique.

TERACEL® Original Haemostat bio-compatible. ISO – 10993 Bio-compatibility results are available to our customers upon request.

TERACEL® Fabric is densely woven knit to give the more powerful ability to absorb heavy blood. Because of similar structure of knit, TERACEL® Fabric is also used as matrix for platelet adhesion and aggregation. The major difference compared to TERACEL® Original is 3 times denser structure and thickness of knit to offer the improved the absorption ability. It is easy-to-use, safe, and effective in controlling bleeding and may help to avoid wound complications as the stemotomy site. In addition, it can be used in wrapping or suturing in place to stop the bleeding and be used in surgery with laparoscopy.

- Excellent Haemostasis in 2 – 3 minutes
- Low pH has antimicrobial effect
- Very good handling characteristics
- Does not stick to instruments
- Typically dissolves in 7 – 14 days
- Bactericidal properties
- Compliance with the requirements of USA
- Plant based product with minimum tissue reaction
- Do not cause irritation to applied body/tissue
- TERACEL® Fabric Haemostat are sterilized by Gamma irradiation technique.

TERACEL® Fabric Haemostat bio-compatible. ISO – 10993 Bio-compatibility results are available to our customers upon request.



PRODUCT	DIMENSIONS	XREF CODE	PACKAGING
TERACEL ORIGINAL	5cm x 35cm	N010535	12 Pcs / Box
TERACEL ORIGINAL	10cm x 20cm	N011020	12 Pcs / Box
TERACEL ORIGINAL	5cm x 7.5 cm	N010575	12 Pcs / Box
TERACEL ORIGINAL	1,25 cm x 5 cm	N011255	12 Pcs / Box
ALCASPONGE ORC	80 mm x 50 mm x 10 mm	HSPONGE	30 Pcs / Box



PRODUCT	DIMENSIONS	X REF CODE	PACKAGING
TERACEL FABRIC	5 cm x 35cm	N020535	12 Pcs / Box
TERACEL FABRIC	10cm x 20cm	N021020	12 Pcs / Box
TERACEL FABRIC	5 cm x 7.5 cm	N020575	12 Pcs / Box
TERACEL FABRIC	1,25 cm x 5 cm	N021255	12 Pcs / Box

TERACEL® Nonwoven is structured of interlocked fibers be needle punching fabric. It is bulky and taken off without difficulty by picking with forceps. It accelerates hemostasis 43% faster by rapid absorbing blood on surgical site in comparison of TERACEL® Original. Non woven structure of TERACEL® provides improved adaptability and adherence in both open and minimally invasive surgery such as cardiovascular, general and peripheral vascular surgery compared to TERACEL® Original. Because it doesn't stick to the sugical instrument, the handling including repositioning and removal during surgery is easily done.

- Excellent Haemostasis in 2 – 3 minutes
- Low pH has antimicrobial effect
- Very good handling characteristics
- Does not stick to instruments
- Typically dissolves in 7 – 14 days
- Bactericidal properities
- Compliance with the requirements of USA
- Plant based product with minimum tissue reaction
- Do not cause irritation to applied body/tissue
- TERACEL® Nonwoven Haemostat are sterilized by Gamma iraddiation technique.

TERACEL® Nonwoven Haemostat bio-compatible. ISO – 10993
Bio-compatibility results are available to our customers upon request.

Teracel® Nonwoven



PRODUCT	DIMENSIONS	XREF CODE	PACKAGING
TERACEL NON-WOWEN	5cm x 35cm	N030535	12 Pcs / Box
TERACEL NON-WOWEN	10cm x 20cm	N031020	12 Pcs / Box
TERACEL NON-WOWEN	5cm x 7,5 cm	N030575	12 Pcs / Box
TERACEL NON-WOWEN	1,25 cm x 5 cm	N031255	12 Pcs / Box

TERACEL® Fibrillar is made from multilayered Oxidized Regenerated Cellulose (ORC), which the surgeon can easily peel with forceps in the desired amount. Each layer in TERACEL® Fibrillar is lightweight and voluminous, designed to melt at the bleeding site and stay in place until the operator removes it. This structure of TERACEL® is suitable for irregular surfaces or tissues, and its voluminous nature is beneficial in orthopedic/spinal and neurological surgery.

Excellent hemostasis within 2-3 minutes.
Antimicrobial effects against various Gram(+), Gram (-), and antibiotic-resistant microorganisms due to its low pH.
Good handling characteristics during surgical operations.
It does not adhere to surgical instruments.
Absorbed within 7-14 days.
Possesses bactericidal properties.
Complies with USP requirements.
Plant-derived product that minimizes tissue reaction.
Does not cause irritation in the applied body/tissue.
TERACEL® Fibrillar Hemostat (Oxidized Regenerated Cellulose) is sterilized using gamma radiation.

TERACEL® Fibrillar Hemostat (Oxidized Regenerated Cellulose) is biocompatible, and ISO 10993 test reports are available.

Teracel® Fibrillar



PRODUCT	DIMENSIONS	REF CODE	PACKAGING
TERACEL® Fibrillar	5cm x 35cm	N040535	12 Pcs / Box
TERACEL® Fibrillar	10cm x 20cm	N041020	12 Pcs / Box
TERACEL® Fibrillar	5cm x 7,5 cm	N040575	12 Pcs / Box
TERACEL® Fibrillar	1,25 cm x 5 cm	N041255	12 Pcs / Box

HEMOSTATIC SPONGE



Hemostatic sponges are effective for both minor and significant bleeding and are often used in situations where traditional methods like sutures or clamps are not feasible. They are eventually absorbed by the body or removed after the bleeding is controlled.



SKIN STAPLER

Skin stapler is a medical device used to close wounds by applying metal staples to the skin. It provides an efficient alternative to traditional suturing, especially in situations where speed is essential. The stapler works by compressing the edges of a wound together and inserting metal staples to hold the tissue in place as it heals. The main advantage of skin staplers is their ability to close wounds quickly and securely, reducing the time the patient spends under anesthesia during surgery or in the emergency room.



BONE WAX

Bone Wax is a sterile, non-absorbable, and waxy substance used primarily in surgical procedures to control bleeding from bone surfaces. It is composed of a blend of beeswax and softening agents, allowing it to be easily molded and applied to bleeding bone edges.





VETERINARIAN SURGERY PRODUCTS

Cassette Sutures • Anti-Bacterial Sutures

CASSETTE SUTURES

Cassette Sutures refer to suture materials that are provided in a spool or cassette format, allowing the veterinarian to cut the desired length as needed during a procedure. Unlike pre-cut and pre-attached needle sutures, cassette sutures offer flexibility and are often used in situations where a variety of suture lengths may be required.



CATGUT CHROME / PLAIN

BRAND	THREAD MATERIAL	EP	USP	LENGTH
WEGUT	Catgut Chrome / Plain	EP3	USP 3/0	100m
WEGUT	Catgut Chrome / Plain	EP4	USP 0	100m
WEGUT	Catgut Chrome / Plain	EP5	USP 1	75m
WEGUT	Catgut Chrome / Plain	EP6	USP 2	50m
WEGUT	Catgut Chrome / Plain	EP7	USP 3	50m
WEGUT	Catgut Chrome / Plain	EP8	USP 4	30m
WEGUT	Catgut Chrome / Plain	EP3.5	USP 2/0	100m

PGA

BRAND	THREAD MATERIAL	EP	USP	LENGTH
WESORB	PGA	EP3	USP 2/0	25m
WESORB	PGA	EP3	USP 2/0	50m
WESORB	PGA	EP5	USP 2	50m
WESORB	PGA	EP5	USP 2	100m
WESORB	PGA	EP6	USP 3	50m
WESORB	PGA	EP6	USP 3	100m
WESORB	PGA	EP7	USP 5	50m
WESORB	PGA	EP3.5	USP 0	25m
WESORB	PGA	EP3.5	USP 0	50m

POLYAMIDE

BRAND	THREAD MATERIAL	EP	USP	LENGTH
WELON	POLYAMIDE	EP3	USP 2/0	100m
WELON	POLYAMIDE	EP6	USP 3	50m
WELON	POLYAMIDE	EP6	USP 3	100m
WELON	POLYAMIDE	EP3.5	USP 0	100m

PDO


BRAND	THREAD MATERIAL	EP	USP	LENGTH
WEDOX	PDO	EP4	USP 1	50m
WEDOX	PDO	EP5	USP 2	50m

PGCL

BRAND	THREAD MATERIAL	EP	USP	LENGTH
WECAPRO	PGCL	EP4	USP 1	50m
WECAPRO	PGCL	EP5	USP 2	50m

ANTIBACTERIAL SUTURES

Antibacterial sutures are surgical sutures that are coated with an antibacterial agent, typically triclosan, to reduce the risk of surgical site infections (SSIs). These sutures are designed to inhibit bacterial colonization on the suture material, which is a potential entry point for infections during the healing process. Antibacterial sutures are particularly beneficial in patients with compromised immune systems or those undergoing surgeries in infection-prone areas.

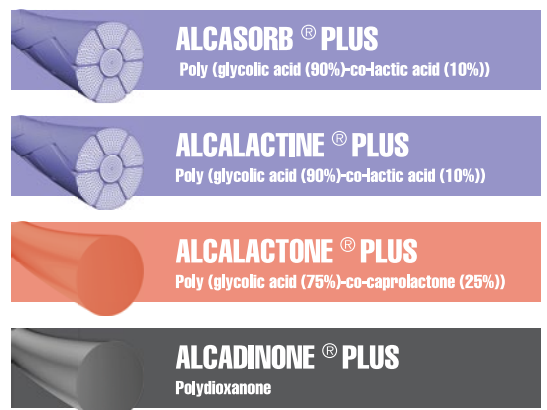



SURGICAL SUTURES

Katsan Surgical Sutures, a leading synthetic surgical suture manufacturing company in Turkey, was founded in 1976 as a catgut producer. Katsan is striving for higher quality products and a bigger market share by using innovative technologies.

With our fully CE marked and ISO 13485 quality certificated product array that are produced in our T.R. Ministry of Health approved facilities, we produce world-class sterile surgical medical devices for a variety of surgical operations.



ABSORBABLE SUTURES

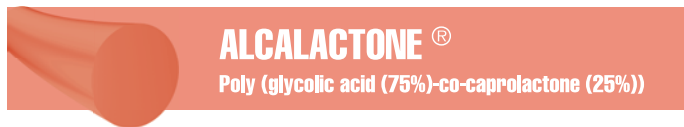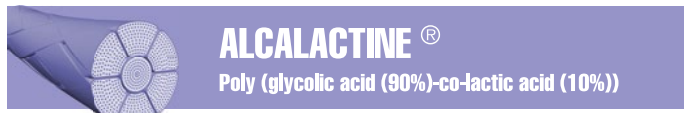

NON-ABSORBABLE SUTURES



AREAS OF USAGE

Plastic surgery, urology, gastrointestinal surgery, skin closure, orthopedic surgery, cholecystectomy, hernia, laparoscopic incisions, peritoneum, fascia, ligation, general tissue approximation, emergency surgery, oral surgery, episiotomy, pediatric surgery, minor surgical procedures, subcutaneous, gastric surgery, intestine, colon, vesica urinaria, vaginal canal, uterus, etc.

PROBIOTICS

* It contains intensively Lactobacillus acidophilus CECT 4529 and Enterococcus Faecium NCIMB 10415 high quality and patented inulin. It also contains Taurine and carnitine and helps heart health.

* Once all the bacteria are encapsulated, it progresses to the large intestines, where it acts, without experiencing any deterioration.

* For the first time in our country, apart from Enterococcus Faecim, a high spectrum of action of Lactobacillus Acidophilus, which has a very important function for dog and cat metabolism, can be provided.

AREAS OF USAGE

Regulating Methabolism

URINAL CURE

POTENT PHYTOTHERAPEUTIC THAT PROMISES TO PROVIDE AN EFFECTIVE ALTERNATIVE TREATMENT METHOD FOR URINARY SYSTEM SYMPTOMS AND DISTAL URETERAL STONES.

AREAS OF USAGE

Urinary system, especially lower urinary tract diseases, constitute a significant problem first in cats and then in dogs. Urolithiasis has a very high incidence, accounting for up to 23% of lower urinary tract diseases in cats and 20.61% in dogs. There are four main types of minerals in cats and dogs; urate, cystine, magnesium ammonium phosphate (struvite) and calcium.

In cats, 68.8% of all uroliths consist of calcium oxalate and 24.2% systruvite, while in dogs 43.6% of all uroliths are struvite and 47% calcium oxalate. Urolithiasis in cats and dogs is a common problem that is still being researched today.

DERMATHOLOGY

Horse chestnut extract, postbiotic, Rice bran oil – gamma oryzanol (sc-co2), Type III collagen, Procapil, Ceramide, D-panthenol, Caffeine + niacinamide, Omega 3-6, Vitamin E, Coconut and palm kernel oil, biotin

AREAS OF USAGE

Procapil suppresses the DHT hormone , which is the main cause of hair loss, and reduces hair loss by ensuring healthy blood flow to the hair follicles and strengthens the hair structure. Provides up to 12% hair increase in 4 weeks of use, anti shedding. Accelerates skin regeneration. Strengthen skin connective tissue. Keeps skin moisture at optimal level. Enriches the lipid wall. Reduces hypersensitivity. Eliminates bad odors. Antibacterial, antifungal, flea repellent, anti-itch. Increases circulation

Antioxidant that prevents dandruff formation. Prevents graying due to melanin deficiency. Strengthens the feather structure. Contains amino acids that enable feather production. Supports keratin formation. Prevents hair from thinning and falling out. Postbiotic content regulates host metabolic pathways and microbiome balance. Strengthens immunity in the skin. Keeps skin oil balanced. Tightens pores. It manages fibroblast levels and increases collagen and elasticity. Protects against the ultraviolet effects of the sun Eliminates hyperpigmentation (after 4 weeks of use). Prevents clumping, conditioner effect

SHAMPOO

Active Ingredients

Horse chestnut extract, tea tree oil, Postbiotic, Rice bran oil- gamma oryzanol (sc-co2), Type III collagen, Procapil, Ceramide, D-panthenol, Caffeine + niacinamide, allantoin, Omega 3-6, Vitamin E, Coconut and palm kernel oil, biotin

AREAS OF USAGE

Procapil suppresses the DHT hormone , which is the main cause of hair loss, and reduces hair loss by ensuring healthy blood flow to the hair follicles and strengthens the hair structure. Provides up to 12% hair increase in 4 weeks of use

anti shedding

Accelerates skin regeneration

Strengthen skin connective tissue

Keeps skin moisture at optimal level

Enriches the lipid wall

Reduces hypersensitivity

Eliminates bad odors

Antibacterial

antifungal

flea repellent

anti-itch

Increases circulation

Antioxidant that prevents dandruff formation

Prevents graying due to melanin deficiency

Strengthens the feather structure

Contains amino acids that enable feather production

Supports keratin formation

Prevents hair from thinning and falling out

Postbiotic content regulates host metabolic pathways and microbiome balance

Strengthens immunity in the skin

Keeps skin oil balanced

Tightens pores

It manages fibroblast levels and increases collagen and elasticity.

Protects against the ultraviolet effects of the sun

Eliminates hyperpigmentation (after 4 weeks of use)

Prevents clumping, conditioner effect

WOUND CARE

Active Ingredients

St. John's wort
Neem oil (rosary tree oil)
Postbiotic lysate
tea tree opioid
Hyaluronic acid with liposome
bovine colostrum
Olive oil
Beeswax
Rice bran oil- gamma oryzanol (sc-co 2)
Ceramide
Type 3 collagen
Lecithin
Tocopherol

It is a natural wound healer known since ancient times.

The polyphenols it contains have antioxidant and anti-inflammatory effects.

The hyperforin enzyme it contains accelerates wound healing by enabling the formation of new blood vessels.

The tannins it contains accelerate tissue regeneration.

Neem oil (rosary tree oil)

It has been used in traditional medicine for wound healing for 2000 years.

It has antibacterial, antifungal and epithelializing effects with the fatty acids, limonoids, vitamin E, triglycerides and antioxidants it contains. Increases collagen production in tissue

Postbiotic lysate

Regulates host metabolic pathways and microbiome balance

Strengthens immunity in the skin

Antimicrobial effect

antioxidant effect

tea tree opioid

Antibacterial, antifungal, antiviral, antiseptic and analgesic effect

Detox effect

anti-inflammatory

Antipruritic effect

Hyaluronic acid with liposome

bovine colostrum

Olive oil

Beeswax

It creates a tissue skeleton and provides the tissue with the energy support necessary for healing.

Significantly shortens the healing process

Rice bran oil- gamma oryzanol (sc-co 2)

Reduces hypersensitivity,

Vitamin E, omega 3-6

Provides anti-inflammatory effect

Increases epithelialization

Increases circulation, accelerates regeneration,

Antioxidant

Ceramide

Ceramides are the main lipid building blocks in the stratum corneum layer of the skin. Provides barrier function of the epidermis

Type 3 collagen

Enables the formation of granulation tissue in the initial phase of wound healing

It provides the necessary support for the increase of new capillaries in the proliferative phase and increases fibroblastic activity at this stage.

Lecithin

Tocopherol

ORAL CARE

Active Ingredients

Tannase (postbiotic) with liposome
Neem oil (rosary tree oil)
Silvia officinalis (sage)
olive oil
Sesame oil
coconut oil
allantoin
Propolis liquid extract
Stevia extract

AREAS OF USAGE

After the spray is applied to the gum, a film quickly forms.

With the liposome technology used, the active ingredients penetrate under the mucosa and act under the tissue for a long time. Thus, maximum long-term therapeutic effect is achieved without loss of active ingredient by swallowing or saliva.

Tannase (postbiotic) with liposome

It increases microcirculation in the gums.
It fights harmful bacteria in the oral flora.
It breaks down dental stones in long-term use.
Prevents bad breath
Neem oil (rosary tree oil)

It is known as the toothbrush tree in India.
It has anti-inflammatory, anti-cavity and anti-tartar properties that accelerate circulation in the gums.
Prevents bad breath
Silvia officinalis (sage)

Preventing pain, inflammation and bad odor, wound healing and antibacterial effect
Sesame oil

Painkiller
whitening effect
coconut oil

Whitens teeth.
It treats and prevents gingivitis and caries.
It kills pathogenic bacteria that cause bad breath.
It treats and prevents bleeding gums.
It strengthens metabolism and provides detoxification.
Reduces tooth sensitivity.
allantoin

keratolytic
Propolis liquid extract

Prevents tooth decay
It disrupts the function of enzymes used by oral bacteria to dissolve sugar.
Helps eliminate tooth sensitivity
Teeth become sensitive due to enamel damage, allowing heat and irritants from the dental canals to enter the dentin, the second layer of the tooth. This causes irritation of the pulp, which is located inside the tooth and where the nerve endings are located. Propolis helps close these channels and eliminates sensitivity. It also supports the regeneration of reparative dentin.
stevia

It is a natural, calorie-free sweetener
Strong antimicrobial effect

oral hygien

BLADDER SUPPORT

Anti-struvite, anti-urolithic, anti inflammatory and anti-bacterial composition.

Complementary against battling calcium oxalate uroliths with its magnesium carbonate content.

Supports regenerating damaged urinary bladder epithelia, assists reducing pain and inflammation.

Helps providing protection against urinary calculus.

Whitens teeth.

It treats and prevents gingivitis and caries.

It kills pathogenic bacteria that cause bad breath.

It treats and prevents bleeding gums.

It strengthens metabolism and provides detoxification.

Reduces tooth sensitivity.

allantoin

keratolytic

Propolis liquid extract

Prevents tooth decay

It disrupts the function of enzymes used by oral bacteria to dissolve sugar.

Helps eliminate tooth sensitivity

Teeth become sensitive due to enamel damage, allowing heat and irritants from the dental canals to enter the dentin, the second layer of the tooth. This causes irritation of the pulp, which is located inside the tooth and where the nerve endings are located. Propolis helps close these channels and eliminates sensitivity. It also supports the regeneration of reparative dentin.

stevia

It is a natural, calorie-free sweetener

Strong antimicrobial effect

oral hygiene

METHABOLIC SUPPORT

LIPOSOMAL GLUTATHION - FELINE

Strong metabolic support with Liposomal Glutathione, Ester C and Vitamin E. Helps with detoxification, regeneration and mitochondrial energy metabolism. Properties that aid the treatment of cancer.

Treatment support for all processes related to oxidative stress. Edible aid in liver diseases, stress; viral, fungal and infectious diseases!

LIPOSOMAL GLUTATHION - CANINE

Functional treat with Liposomal Glutathione, Ester C, Vitamin E and Alpha Lipoic Acid(ALA). ALA is a powerful antioxidant that removes toxic substances and supports the treatment of diabetes! Properties that aid the treatment of cancer.

Treatment support for all processes related to oxidative stress. Edible aid in the shape of treats that support the treatment of liver diseases, stress; viral, fungal and infectious diseases!

HAIRBALL TREATMENT

HAIRBAL CURE

Edible intestinal satiety product powered by postbiotic lysate and psyllium.

Support for balancing the gut microbiome, easing digestion and also regulating numerous immune and metabolic functions. Helps optimal microbiome composition, enhanced immune function, hormonal pathways, metabolic functions, and even general mood.